

PAA310Hu01

Polyclonal Antibody to Apolipoprotein H (APOH)
Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: APOH, Human **Purification:** Affinity Chromatography.

Clonality: Polyclonal Applications: WB, ICC, IHC-P, IHC-F, ELISA

Host: Rabbit Concentration: 200µg/mL

Immunoglobulin Type: IgG **UOM**: 100μg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant APOH (Thr22~Cys345) expressed in *E.coli*.

Accession No.: RPA310Hu01

Sequence: The target protein is fused with N-terminal His-Tag and its sequence

is listed below.

MGHHHHHHSGSEF-TCPKPDDLP FSTVVPLKTF YEPGEEITYS CKPGYVSRGG MRKFICPLTG LWPINTLKCT PRVCPFAGIL ENGAVRYTTF EYPNTISFSC NTGFYLNGAD SAKCTEEGKW SPELPVCAPI ICPPPSIPTF ATLRVYKPSA GNNSLYRDTA VFECLPQHAM FGNDTITCTT HGNWTKLPEC REVKCPFPSR PDNGFVNYPA KPTLYYKDKA TFGCHDGYSL DGPEEIECTK LGNWSAMPSC KASCKVPVKK ATVVYQGERV KIQEKFKNGM LHGDKVSFFC KNKEKKCSYT EDAQCIDGTI EVPKCFKEHS SLAFWKTDAS DVKPC

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against APOH. It has been selected for its ability to recognize APOH in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500




Immunohistochemistry in formalin fixed frozen section: 1:100-500


Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[IMAGES]

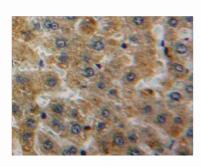

Used in DAB staining on fromalin fixed paraffin-embedded Liver tissue

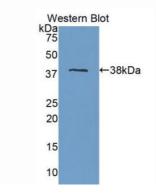

Used in Western Blot ,Sample:
Recombinant APOH.Human

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant APOH (Thr22~Cys345) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN₃ 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.