

#### PAA541Po01

Polyclonal Antibody to Brain Natriuretic Peptide (BNP)

Organism Species: Sus scrofa; Porcine (Pig)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

## [ PRODUCT INFORMATION ]

Immunogen: BNP, Porcine

Clonality: Polyclonal

Host: Rabbit


Immunoglobulin Type: IgG

Purification: Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

**UOM**: 50μg

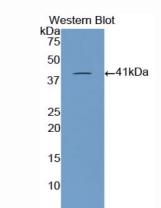

Sample: Recombinant BNP, Porcine

#### [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant BNP (His26~Tyr131) expressed in E.coli.

Accession No.: RPA541Po01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and

GST-tag and its sequence is listed below.

MSPILGYWKI KGLVQPTRLL LEYLEEKYEE HLYERDEGDK WRNKKFELGL EFPNLPYYID GDVKLTQSMA IIRYIADKHN MLGGCPKERA EISMLEGAVL DIRYGVSRIA YSKDFETLKV DFLSKLPEML KMFEDRLCHK TYLNGDHVTH PDFMLYDALD VVLYMDPMCL DAFPKLVCFK KRIEAIPQID KYLKSSKYIA WPLQGWQATF GGGDHPPKSD GSTSGSGHHH HHHSAGLVPR GSTAIGMKET AAAKFERQHM DSPDLGTLEV LFQGPLGSEF- HPLGG AGLASELPGI QELLDRLRDR VSELQAERTD LEPLRQDRGL TEAWEAREAA PTGVLGPRSS IFQVLRGIRS PKTMRDSGCF GRRLDRIGSL SGLGCNVLRR Y



## [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against BNP. It has been selected for its ability to recognize BNP in immunohistochemical staining and western blotting.

### [APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

## [CONTENTS]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

# [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant BNP (His26~Tyr131) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate. 5uL per well when used in enhanced chemilumescent (ECL).

**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.

### [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.