### PAB466Hu01 Polyclonal Antibody to Connexin 37 (CX37) **Organism Species: Homo sapiens (Human)** Instruction manual


#### FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES


[ PRODUCT INFORMATION ] Immunogen: CX37, Human **Clonality:** Polyclonal Host: Rabbit Immunoglobulin Type: IgG

9th Edition (Revised in Jul, 2013)

Sample: Recombinant CX37, Human

Sequence: The target protein is fused with N-terminal His-Tag and its sequence MGHHHHHHSG SEF- VL QFLFVSTPTL VYLGHVIYLS RREERLRQKE GELRALPAKD PQVERALAAV ERQMAKISVA EDGRLRIRGA LMGTYVASVL CKSVLEAGFL YGQWRLYGWT MEPVFVCQRA PCPYLVDCFV SRPTEKTIFI IFMLVVGLIS LVLNLLELVH LLCRCLSRGM RARQGQDAPP TQGTSSDPYT DQVFFYLPVG QGPSSPPCPT YNGLSSSEQN WANLTTEERLASSRPPLFLD PPPQNGQKPP SRPSSSAS

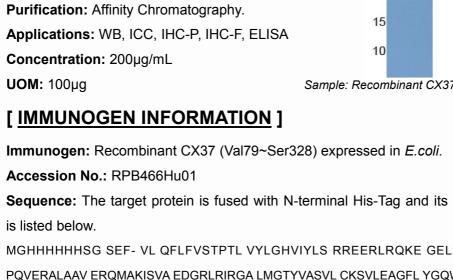

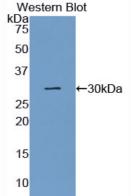

# [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against CX37. It has been selected for its ability to recognize CX37 in immunohistochemical staining and western blotting.

### [APPLICATIONS]

Western blotting: 1:100-400 Immunocytochemistry in formalin fixed cells: 1:100-500 Immunohistochemistry in formalin fixed frozen section: 1:100-500 Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

# [<u>CONTENTS</u>]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

## [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant CX37 (Val79~Ser328) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$  0.02%.

# [<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.