
#### PAB932Mu01 Polyclonal Antibody to Glyceraldehyde-3-Phosphate Dehydrogenase (GAPDH) Organism Species: Mus musculus (Mouse) Instruction manual

#### FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

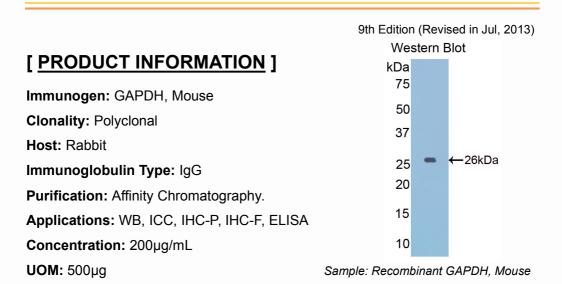

### [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant GAPDH (Gly56~Phe284) expressed in E.coli.

Accession No.: RPB932Mu01

**Sequence:** The target protein is fused with N-terminal His-Tag and its sequence is listed below.

MGHHHHHHSGSEF- GTVKA ENGKLVINGK PITIFQERDP TNIKWGEAGA EYVVESTGVF TTMEKAGAHL KGGAKRVIIS APSADAPMFV MGVNHEKYDN SLKIVSNASC TTNCLAPLAK VIHDNFGIVE GLMTTVHAIT ATQKTVDGPS GKLWRDGRGA AQNIIPASTG AAKAVGKVIP ELNGKLTGMA FRVPTPNVSV VDLTCRLEKP AKYDDIKKVV KQASEGPLKG ILGYTEDQVV SCDF

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008 ; ISO13485:2003 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.doud-clone.us | E-mail: mail@cloud-clone.us
Export Processing Zone Building F, Wuhan, Hubel 430056, PRC | Toll free: 0086-800-880-0667 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@cloud-clone.us

# [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against GAPDH. It has been selected for its ability to recognize GAPDH in immunohistochemical staining and western blotting.

### [APPLICATIONS]

Western blotting: 1:50-400 Immunocytochemistry in formalin fixed cells: 1:50-500 Immunohistochemistry in formalin fixed frozen section: 1:50-500 Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

# [<u>CONTENTS</u>]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

## [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant GAPDH (Gly56~Phe284) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$  0.02%.

# [<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.cloud-clone.us | E-mail: mail@cloud-clone.us Export Processing Zone Building F, Wuhan, Hubei 430056, PRC | Toll free: 0086-800-880-0667 | Fax: 0066-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@cloud-clone.us