
#### PAD982Hu01 Polyclonal Antibody to Dimethylarginine Dimethylaminohydrolase 1 (DDAH1) Organism Species: Homo sapiens (Human) Instruction manual

#### FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

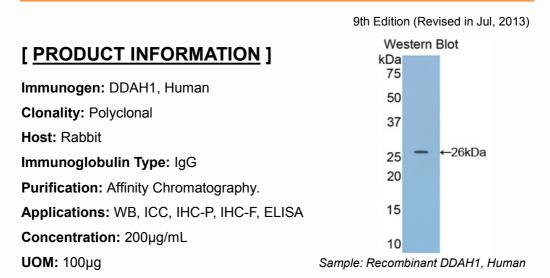

#### [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant DDAH1 (Ala20~Pro215) expressed in *E.coli*. Accession No.: RPD982Hu01 Sequence: The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below. MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF-A LPESLGQHAL RSAKGEEVDV ARAERQHQLY VGVLGSKLGL QVVELPADES LPDCVFVEDV AVVCEETALI TRPGAPSRRK EVDMMKEALE KLQLNIVEMK DENATLDGGD VLFTGREFFV GLSKRTNQRG AEILADTFKD YAVSTVPVAD GLHLKSFCSM AGPNLIAIGS SESAQKALKI MQQMSDHRYD KLTVP

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008; ISO13485:2003 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toil free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.doud-clone.us | E-mail: mail@cloud-clone.us Export Processing Zone Building F, Wuhan, Hubel 430056, PRC | Toil free: 0086-800-880-0687 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@uscnk.com

## [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against DDAH1. It has been selected for its ability to recognize DDAH1 in immunohistochemical staining and western blotting.

### [APPLICATIONS]

Western blotting: 1:50-400 Immunocytochemistry in formalin fixed cells: 1:50-500 Immunohistochemistry in formalin fixed frozen section: 1:50-500 Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

# [<u>CONTENTS</u>]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

## [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant DDAH1 (Ala20~Pro215) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$  0.02%.

# [<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.