
PAB255Mu01 Polyclonal Antibody to Elongin A (ELOA) Organism Species: Mus musculus (Mouse) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

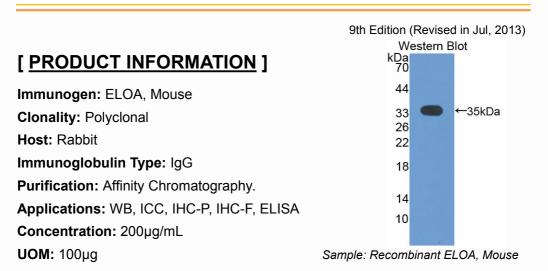

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant ELOA (Val76~Ser345) expressed in E.coli.

Accession No.: RPB255Mu01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- VPVER NSEAEDQDFE KNNSRKRPRD ALQREEELEG NYQESWKPSG SRSYSPEHRQ KKHKKLSEPE RPHKVAHSHE KRDERKRCHK VSPPYSSDPE SSDYGHVQSP PPSSPHQMYT DLSRSPEEDQ EPIISHQKPG KVHSNTFQDR LGVSHLGEQG KGAVSHHKQH RSSHKEKHPA DAREDEKISA VSREKSHKAS SKEESRRLLS GDSAKEKLPS SVVKKDKDRE GSSLKKKFSP ALDVASDNHF KKPKHKDSEK AKSDKNKQSV DGVDS

Designed by Cloud-Clone Corp., Assembled by Uscn Life Science Inc. ISO9001:2008 ; ISO13485:2003 11271 Richmond Avenue Suite H104, Houston, TX 77082, USA | Toll free: 001-888-960-7402 | Fax: 001-832-538-0088 | Http://www.doud-clone.us | E-mail: mail@cloud-clone.us
Export Processing Zone Building F, Wuhan, Hubel 430056, PRC | Toll free: 0086-800-880-0667 | Fax: 0086-27-8425-9551 | Http://www.uscnk.com | E-mail: mail@uscnk.com

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against ELOA. It has been selected for its ability to recognize ELOA in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400 Immunocytochemistry in formalin fixed cells: 1:100-500 Immunohistochemistry in formalin fixed frozen section: 1:100-500 Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

[<u>CONTENTS</u>]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant ELOA (Val76~Ser345) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$ 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.