PAA468Hu01 Polyclonal Antibody to Decay Accelerating Factor (DAF) Organism Species: Homo sapiens (Human) *Instruction manual*

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: DAF, Human Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant DAF (Asp35~Gly285) expressed in E.coli.

Accession No.: RPA468Hu01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below.

MGSSHHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- DCGLPP DVPNAQPALE GRTSFPEDTV ITYKCEESFV KIPGEKDSVI CLKGSQWSDI EEFCNRSCEV PTRLNSASLK QPYITQNYFP VGTVVEYECR PGYRREPSLS PKLTCLQNLK WSTAVEFCKK KSCPNPGEIR NGQIDVPGGI LFGATISFSC NTGYKLFGST SSFCLISGSS VQWSDPLPEC REIYCPAPPQ IDNGIIQGER DHYGYRQSVT YACNKGFTMI GEHSIYCTVN NDEGEWSGPP PECRG

[ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against DAF. It has been selected for its ability to recognize DAF in immunohistochemical staining and western blotting.


[APPLICATIONS]

Western blotting: 1:100-400 Immunocytochemistry in formalin fixed cells: 1:100-500 Immunohistochemistry in formalin fixed frozen section: 1:100-500

Cloud-Clone Corp.

Immunohistochemistry in paraffin section: 1:50-200 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.


[<u>IMAGES</u>]

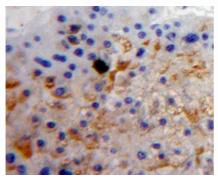

Used in DAB staining on fromalin fixed

paraffin-embedded Liver tissue

[CONTENTS]

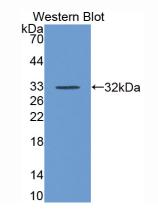

Used in Western Blot, Sample: Recombinant DAF, Human

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant DAF (Asp35~Gly285) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN $_3$ 0.02%.

[<u>STORAGE</u>]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.